


# Ashley and Bovingdon Reflex Balloon Catheters

Practitioners in anorectal physiology now have a choice of two catheters custom-designed for balloon distension, R.A.I.R., expulsion and other common tests.

## *The Ashley reflex balloon catheter (ARC 2410/10 & ARC 2410/100)*

Has a 400mm PVC tube-section securely fitted with a large, thin-walled latex balloon for maximum physiological compliance, which may be inflated up to 1,000ml if needed. A 3-way tap is fitted to the proximal end of the catheter, to enable easy and controlled inflation/deflation in stages.

The catheter is available in a box of 10, individually peel-pouched, or in a box of 100 bulk-packed for easier handling in a busy clinic.

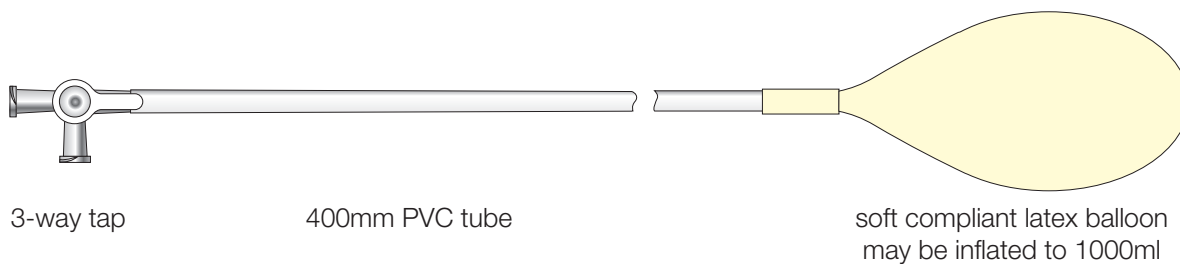


## *The Bovingdon silicone balloon catheter (ARC 2470 & ARC 2475)*

In addition to distension, RAIR and expulsion testing, is also suitable for both intra-vaginal and anorectal pressure biofeedback training.

The 70mm shaped silicone balloon is fitted to either a 400mm or 1000mm PVC catheter tube, and a proximal 3-way tap. The balloon may be inflated with air or water up to a capacity of 500ml.

The Bovingdon ARC catheter is available in a box of 10, individually peel-pouched, or in a box of 25x catheter packs each containing a 1000mm catheter and 50ml syringe.

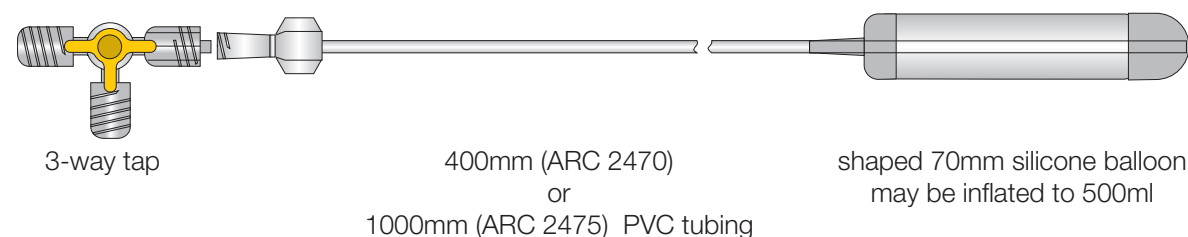

# Ashley and Bovingdon Reflex Balloon Catheters



## *The Ashley*

- 400mm-long PVC catheter tubing
- Large 1,000ml capacity bendy and physiologically compliant balloon
- Designed for RAIR, Distension, Expulsion & other tests
- Balloon is low protein yield latex
- Balloon is white, to assist in the visual diagnosis of pathologies
- 3-way tap fitted as standard to all catheters
- Industrially manufactured to comply with ISO13485:2003, MDD93/42/EEC, CE-Marked
- Single use, non-sterile

### Available in

- Box of 10, individually peel-pouched: order code ARC 2410/10
- Box of 100, bulk-packed for ease of use in busy clinic: order code ARC 2410/100

## *The Bovingdon*

- 400mm or 1000mm-long PVC catheter tubing
- Shaped, cylindrical, silicone balloon
- Latex-free
- Balloon is 70mm long, with 500ml fill-capacity
- Use for RAIR, Distension, Expulsion & other GI Physiology tests
- Suitable for intra-vaginal or ano-rectal biofeedback training
- Industrially manufactured to comply with ISO13485:2003, MDD93/42/EEC, CE-Marked
- Single use, non-sterile

### Available in

- Box of 20, 400mm catheter variant, individually peel-pouched: order code ARC 2470
- Box of 25, 1000mm catheter variant in procedure pack, with 3-way tap and 50ml syringe: order code ARC 2475


order code	description	catheters per box
ARC 2410/10	Ashley reflex balloon catheter: 400mm catheter tubing with 1 litre soft latex balloon & 3-way tap. Individually peel-pouched	10
ARC 2410/100	Ashley reflex balloon catheter: 400mm catheter tubing with 1 litre soft latex balloon & 3-way tap. Double bagged in bulk for easy handling in clinic	100
ARC 2470	Bovingdon reflex balloon catheter: 400mm catheter tubing with 70mm silicone balloon & 3-way tap. Individually peel-pouched	20
ARC 2475	Bovingdon reflex balloon catheter: 1000mm catheter tubing with 70mm silicone balloon. Packed with 50ml syringe & 3-way tap	25



Please contact us for information on the complete catheter range, for samples, or to order.

Tel: +44 (0)1494 794315  
Fax: +44 (0)1494 772759  
info@stericom.com  
www.stericom.com

Stericom Ltd  
Units 1&2 Higham Mead  
Chesham HP5 2AH England

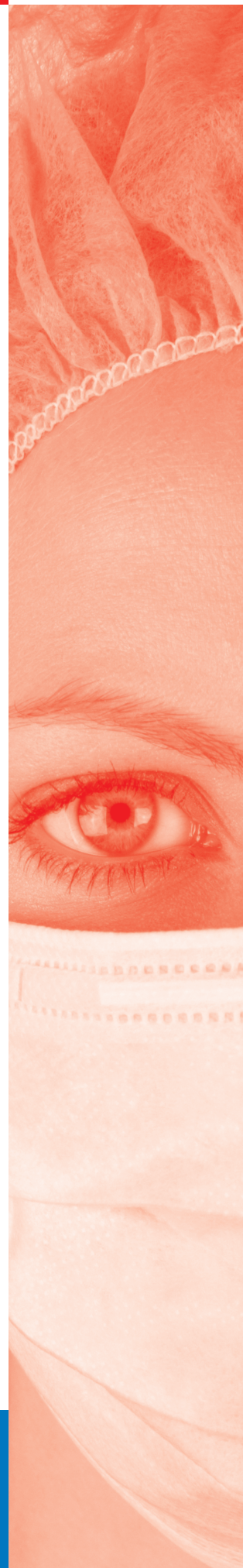


# Langley Diagnostic and Bolton Trainer Catheters

Physiotherapists and other practitioners specializing in women's health now have a range of catheters designed for the effective diagnosis and treatment of fecal incontinence. The 2-balloon Langley Catheter may be used for RAIR, resting/squeeze pressures and associated tests without resorting to multi-channel manometry. The Bolton trainer catheters are custom-designed for easy and effective pressure biofeedback therapy.

## *The Langley 2-lumen, twin-balloon diagnostic catheter (ARC 2200)*


This diagnostic catheter has two balloons fixed in tandem on one 2-channel shaft, each connected to a 1500mm color-coded pressure line. Each line is Ch6 and fitted with a 3-way tap to enable rapid inflation and deflation.

A number of investigations may be performed easily by the practitioner with this catheter: sphincter resting and squeeze pressures, simplified RAIR test without the need for multi-channel water-perfused manometry, measurement of sensitivity threshold and monitoring progress whilst teaching sphincter re-education exercises. Full protocols are supplied in the instructions for use in each box.

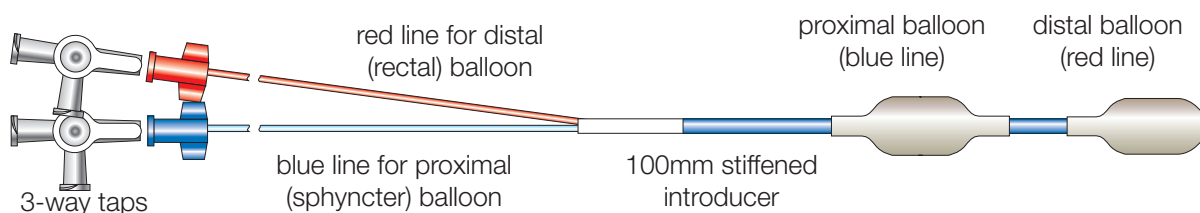

## *The Bolton biofeedback trainer catheter (ARC 2100 & ARC 2101)*

Designed for anorectal pressure biofeedback retraining, the Bolton Catheter has 2000mm Ch6 PVC tubing to enable rapid inflation, deflation and response, and is fitted with a luer connector and 3-way tap. Either a green PVC or latex balloon is available and both are fitted with a 100mm stiffened introducer to enable easy insertion by either the practitioner or the patient.

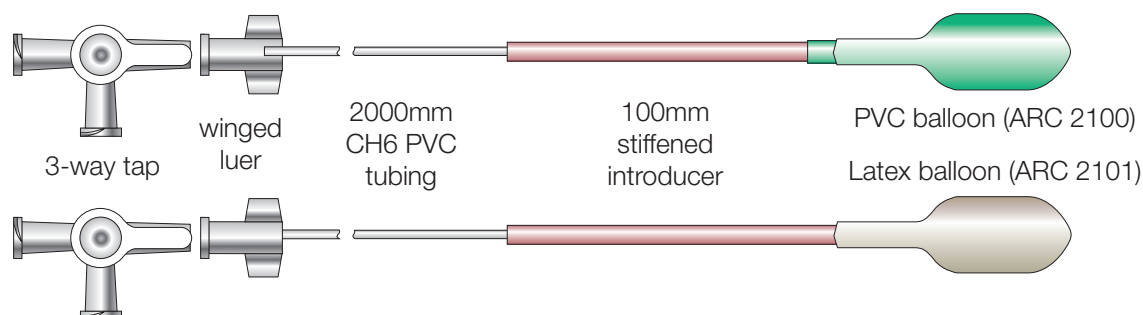

# Langley Diagnostic and Bolton Trainer Catheters



## The Langley

- Two latex balloons fitted in tandem on a Ch14, 2-channel shaft
- Shaft is stiffened but also bendy and compliant to enable easy insertion with maximum patient comfort
- 2x color-coded Ch6, 1500mm long pressure lines
- Each pressure line fitted with 3-way tap
- Use for resting/squeeze pressures, RAIR, measuring sensitivity to distension, expulsion testing and other investigations
- Industrially manufactured to comply with ISO13485:2003, MDD93/42/EEC, CE-marked
- Single use, non-sterile

### Available in

- Box of 5, individually peel-pouched: order code ARC 2200. Instructions for use/investigation protocols included in each box

## The Bolton

- 2000mm Ch6 catheter tubing for ease of inflation/deflation
- Luer connector and 3-way tap
- Either PVC or latex balloon for sphincter retraining
- 100mm stiffened introducer incorporated into catheter/balloon intersection to enable easy placement by practitioner or patient
- Use for pressure biofeedback retraining
- Industrially manufactured to comply with ISO13485:2003, MDD93/42/EEC, CE-marked
- Single use, non-sterile

### Available in

- Box of 20, individually peel-pouched: order code ARC 2100 with green PVC balloon
- Box of 20, individually peel-pouched: order code ARC 2101 with latex balloon

order code	description	catheters per box
ARC 2100	Bolton biofeedback trainer catheter. 2000mm with PVC balloon, introducer & 3-way tap	20
ARC 2101	Bolton biofeedback trainer catheter. 2000mm with latex balloon, introducer & 3-way tap	20
ARC 2200	Langley 2-lumen, twin-balloon anorectal catheter, includes 3-way tap on each line	5



Please contact us for information on the complete catheter range, for samples, or to order.

Tel: +44 (0)1494 794315  
Fax: +44 (0)1494 772759  
info@stericom.com  
www.stericom.com

Stericom Ltd  
Units 1&2 Higham Mead  
Chesham HP5 2AH England

